



SUNDAL IMPEX

LIPOSUCTION CANNULAS & ACCESSORIES



Sundal Impex is a manufacturer and supplier of fine surgical instruments. It has especially been designed for the purpose of supplying and providing the medical profession with high quality and state of the art products/ instruments. It is continuously developing new products, upgrading its instruments in all possible ways through technological enhancement and investigating under improved methods.

Sundal Impex welcomes and appreciates queries and inquiries regarding instruments whether shown in our catalog or not. We are convinced that the patterns shown are widely used and accepted. A variant style is also preferable according to demand and feasibility. In such cases, we can supply the instrument on catalog number, sample or a brief description is provided.

Table Of Content

INFILTRATION CANNULAS	1
GENERAL SUCTION CANNULAS	2
HARVESTING CANNULAS	3
SPECIALITY CANNULAS	4
CANNULA HUBS / CONNECTIONS	5-7

Tumescent Infiltrator**01**

This cannula is used specifically for the infusion of tumescent solution. The configuration of the ports allows even distribution of the fluid.

**Tumescent Infiltrator Angled****02**

This cannula is used specifically for the infusion of tumescent solution. The configuration of the ports allows even distribution of the fluid.

**Facelift Infiltration Cannula****03**

This cannula is used specifically for the infusion of tumescent solution. The configuration of the ports allows even distribution of the fluid.

**GENERAL SUCTION CANNULAS****Stevens Style Cannula****04****Cobra Bibevel Style****05**

More aggressive than the Standard, Pyramid or Mercedes, this cannula has a tapered tip with two ports and works

**Cobra Style Round Tip****06**

Less aggressive than the Cobra Bibevel this cell friendly cannula has a rounded tip with two ports and works effectively within dense fibrous tissue during liposuction procedures.



Standard, Single Hole
07

This is the most common cell friendly cannula design for soft tissue removal. The single port is located 180 degrees from the thumb rest on the hub.


Pyramid Holes Cannula
08

This commonly used design has three ports located 180 degrees from the thumb print. It is a slightly more aggressive cannula than the Standard during liposuction procedures.


Mercedes
09

This circumferential three port configuration cell friendly liposuction cannula is effective for multi-plane soft tissue harvesting.


3 Spiral Port Cannula (Carraway Style)
10

Spiral three port design on this Cell Friendly liposuction cannula is effective for multi-plane adipose tissue harvesting.


Six Spiral Port (Adjust below 3 Spiral Ports)
11

Six port spiral configuration on this Cell Friendly liposuction cannula harvests maximum tissue volume in minimal time.


Three Port Cannula
12

This body sculpting instrument is excellent for harvesting small aliquots of tissue during liposuction procedures.


Basket Shaped Cannula
13

The Cell Friendly Basket Cannula is a versatile liposculpting tool. It can operate as an equalizer without vacuum to treat adipose tissue and skin irregularities in revision liposuction and fat transfer procedures. The Basket Cannula is also adept under vacuum as a fat harvesting cannula. For finishing, this cannula can be used as an injector for equalizing and adding adipose tissue to sculpt smoother contours.



Multiport Speed Harvester- Both End Sharp
14

A multiport harvester with elevated 1mm microports which yield small enough aliquots of fat tissue to allow for easy and smooth injection with the Tulip 0.7, 0.9 and 1.2mm injectors.


Multiport Speed Harvester- One End Sharp
15

A multiport harvester with elevated 1mm microports which yield small enough aliquots of fat tissue to allow for easy and smooth injection with the Tulip 0.7, 0.9 and 1.2mm injectors.


Multi Port Harvester - Both End Bevel
16

A sophisticated, precision-designed harvester that is truly Patient Friendly. The ports in the Sorensen Harvester are machined with alternating 60 degree cutting edges. This unique design is less traumatic than other multi-port harvesters and delivers a fine-grained, high fat content/purity yield that can easily pass through 1.2mm, 0.9mm and 0.7mm blunt injectors.


HARVESTING CANNULAS
17

Fig. 2B
18

Fig. 1C
19

**Coleman Liposuction
Cannula Fig. MIC I**
20

**Coleman Liposuction
Cannula Spoon Tip Fig.
MIC III**

This tip works well for needle extraction of soft tissue and as a soft tissue re injection needle.

21

Little Injector

Open and more pointed blunt tip is very effective in dense and fibrous tissue.

V-Dissector**22**

This cannula is specially designed for use in the superficial compartment of the truncal area. The cutting notch is very useful for undermining in facial procedures and in dense, fibrous tissue.

**Flap****23**

The Flap is primarily used as an undermining dissector, and in smaller diameters, for facial tissue removal. The Flap Cannula in a 2.1mm diameter is very useful for undermining the nasolabial fold prior to autologous tissue injections, as well as for tissue removal in this anatomic area.

**Field Flap****24**

With its elevated leading edge, this cannula is designed to open and expose soft tissue for more precise suctioning along the mandible and sub-mental areas.





LL

Luer lock




UL

Ultra Luer lock



TH

60cc Toomey hub

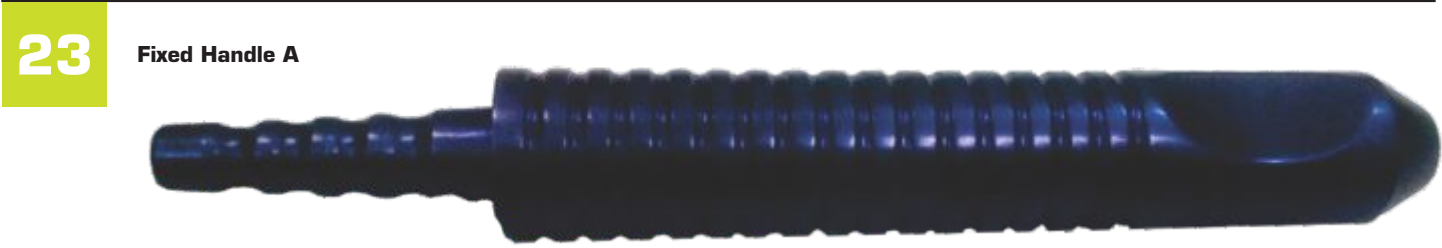

23

Fixed Handle A





24

Fixed Handle B

Sizes Chart

000	Ø 1.0 mm X 10 cm
001	Ø 1.5 mm X 12 cm
002	Ø 2 mm X 15 cm
003	Ø 2 mm X 20 cm
004	Ø 2 mm X 25 cm
005	Ø 2 mm X 30 cm
006	Ø 2 mm X 35 cm
007	Ø 3 mm X 15 cm
008	Ø 3 mm X 20 cm
009	Ø 3 mm X 25 cm
010	Ø 3 mm X 30 cm
011	Ø 3 mm X 35 cm
012	Ø 4 mm X 15 cm
013	Ø 4 mm X 20 cm
014	Ø 4 mm X 25 cm
015	Ø 4 mm X 30 cm
016	Ø 4 mm X 35 cm
017	Ø 5 mm X 15 cm
018	Ø 5 mm X 20 cm
019	Ø 5 mm X 25 cm
020	Ø 5 mm X 30 cm
021	Ø 5 mm X 35 cm
022	Ø 6 mm X 15 cm
023	Ø 6 mm X 20 cm
024	Ø 6 mm X 25 cm
025	Ø 6 mm X 30 cm
026	Ø 6 mm X 35 cm

Order reference procedure

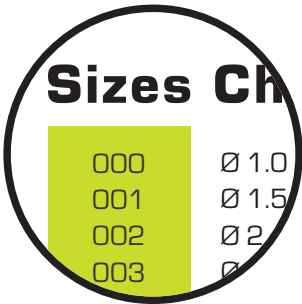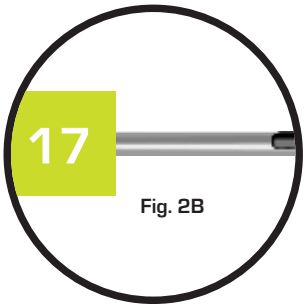

Product Code + HUB Code + Size Chart
Example: 16-LL-000

Power Injector

L-PI-001	20 cc
L-PI-002	10 cc
L-PI-003	Universal (for both 20 and 10cc syringes)

The Power Injector allows the surgeon to re-inject soft tissue through large or small diameter injectors with very little effort and maximum control. The Universal Power Injector will hold a 10cc or a 20cc BD syringe or a BD syringe or 60cc Monoject syringe.


Fat Syringe

L-FS-001	50/60 cc
L-FS-002	20 cc
L-FS-003	10 cc
L-FS-004	5 cc
L-FS-005	3 cc

The fat syringe is an ideal tool to help surgeons for manual liposuction by creating negative pressure.


SnapLoks

L-SL-001	60 cc
L-SL-002	50 cc
L-SL-003	35 cc
L-SL-004	20 cc
L-SL-005	10 cc

These syringe snap locks fit into the lunger angles of the syringe. When the plunger is pulled out the lock snaps open on the barrel lip thus holding the syringe under vacuum.




Needle Connector**L-NC-001** Length 2.3"

Needle connector for use with reusable tubing.

**IV Spike****L-IV-001** Stainless Steel

IV connector for use with reusable tubing.

**Luer to Luer****L-LL-001****Decant Stands****L-DS-001** 60cc to Luer
L-DS-002 Luer to Luer

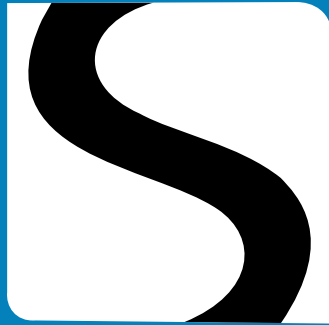

SUNDAL IMPEX

LIPOSUCTION CANNULAS & ACCESSORIES



Roras Road, Muzaffar pur, Sialkot 51310 Pakistan.

@ info@sundalimpex.com

🌐 www.sundalimpex.com

☎ +92523572564